## **HEADQUARTERS**

Apex Medical Corp. No. 9, Min Sheng St., Tu-Cheng, New Taipei City, 23679, Taiwan T +886 2 2268 5568 F +886 2 2268 6525 marketing@apexmedicalcorp.com

## PAN-ASIA

Apex Medical Corp. No. 9, Min Sheng St., Tu-Cheng, New Taipei City, 23679, Taiwan T +886 2 2268 5600 /0800 054 133 F +886 2 2268 6526 customer@apexmedicalcorp.com

### STURDY INDUSTRIAL

168, Sec.1, Zhongxing Road, Wugu District, New Taipei City 24872, Taiwan T +886 2 2981 5245 F +886 2 2984 8208 foreign@sturdy.com.tw www.sturdy.com.tw

### CHINA

Apex Medical (Kunshan) Corp. No. 1368, ZiZhu Road, Kunshan, Jiang Su Province, China T +86 512 8617 1660 Ext 1150 F +86 512 8617 8663 customer.cn@apexmedicalcorp.com

## INDIA

Apex Medicalcorp India Pvt. Ltd. 303, Third Floor Block 4A DLF Corporate Park Phase-ill, Gurgaon, Haryana T +91 11 6654 4138 F+911166544052 enquiries.india@apexmedicalcorp.com

## THAILAND

Apex Medical Thailand 65/100 Chamnanphenjati Business Center Bldg. 11/F Rama 9 Road, Huaykwang Bangkok 10320 T +662 948 4848 F +662 948 4170 sales.thai@apexmedicalcorp.com

## MIDDLE EAST • AFRICA

29, Misr Liltameer Building, Stree#10, AREA#8, Sheraton Building Heliopolis Cairo, Egypt. T +202 2 2671303 F+20222678649 raafat.wahby@apexmedicalcorp.com

## PAN-EUROPE

Apex Medical S.L. Elcano 9, 6ª planta 48008 Bilbao. Vizcaya. Spain T +34 94 4706408 F +34 94 4706409 info@apexmedical.es

Apex Medical Limited Unit 33, Great Western Business Park McKenzie Way, Worcester ,WR4 9GN United Kingdom T +44 (0) 190 5774 695 sales@apexmedicalcorp.co.uk

## FRANCE

Apex Medical France 4 Boulevard de la Chanterie 49124 Saint-Barthelemy d'Anjou France T +33 2 72 79 27 13 F +33 2 85 73 07 67 info.france@apexmedicalcorp.com

SLK Vertriebsgesellschaft GmbH Lindenhorster Straße 38-40 44147 Dortmund T +49 (0)231-925360-175 F +49 (0)231-925360-29 info@slk-gmbh.de www.slk-gmbh.de

## **AMERICAS**

Apex Medical USA 927 Mariner Street, Brea, CA 92821, USA T +1 714 671 3818 F +1 714 494 8173 sales.usa@apexmedicalcorp.com

Distributed by

Print-2019/08-500 All rights reserved



# Domus 3D

## Alternating Pressure Redistribution System

- Suitable for patients in medium to high risk of pressure injury
- Digital pump for precise pressure control
- Key functions to satisfy pressure ulcer therapy and overall prevention
- Continuous Low Pressure Therapy, increased comfort for better care












Apex Medical Corp. Domus 3D

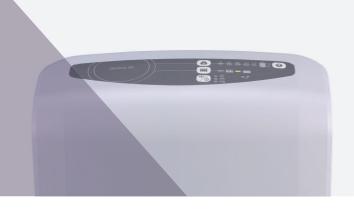


# Domus 3D

Alternating Pressure Redistribution System

"Your essential solution"

## Pressure Relief During Alternating

Domus 3D is an essential solution which not only offers alternating pressure redistribution, but also has a static mode with continuous low pressure

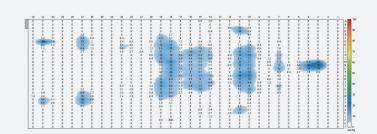




Cell-on-cell design (8") The bottom cells are capable of providing adequate support for the system without an extra foam base in the case of power failure



Low Air Loss Provides good ventilation and reduces the accumulation of heat and moisture



Quick Connector with Cap Provides a stable surface for patient when transporting

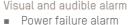

- Low pressure alarm
- Tech. Support alarm

Panel lock



Helps prevent unauthorized changes to the pump settings





Facile foam compartment (5") Extra foam interlayer with buttoned design as an adequate support base



Storage strap Side buckle can simply wrap and store mattress



Nylon/PU Cover Vapour permeable cover offers you the comfort and protection you need









Cable management Storage strap on the side of mattress for managing cable





## Digital for precision

■ Digital pump makes precise pressure control possible for better patient care. Enables the compressor to work more efficiently, creating a longer-lasting system

## Alternating

■ Continuously and sequentially inflate and deflate air cells to mimic natural sleep movement patterns and avoid long term pressurization of tissue, while stimulating blood circulation

## Continuous Low Pressure (CLP)

Reduces the pressure precisely by only 1/3 from the alternate pressure level, providing a static comfort surface while still retaining sufficient support

## Specifications

Model		Domus 3D	
Pump	Dimension	29x18.5x12.6 cm	
	Weight	2.3Kg	
	Case Material	Flame retardant ABS	
	Supply Voltage	AC 220 - 240V / 50Hz ; AC 120V / 60Hz	
Mattress		Overlay	Replacement
	Dimension	190 x 85 x 4"	200 x 80/85/90 x 8"
		200 x 80/85/90 x 5"+3	
	Cell Height	17 x 4"/5" cells	20 x 8" cells (cell-on-cell)
	Weight	5.6 kg	6.9 kg
	Cover Material	Nylon PU (Quilted)	Nylon PU (Quilted)
		4 way Stretch (optional)	4 way Stretch (optional)
	Cell Material	Nylon PU	Nylon PU
	Maximum Patient Weight	200 kg	250 kg
	Flame Retardant Standards	EN597-1; EN597-2	EN597-1; EN597-2

Specifications are subject to change for improvement without notice